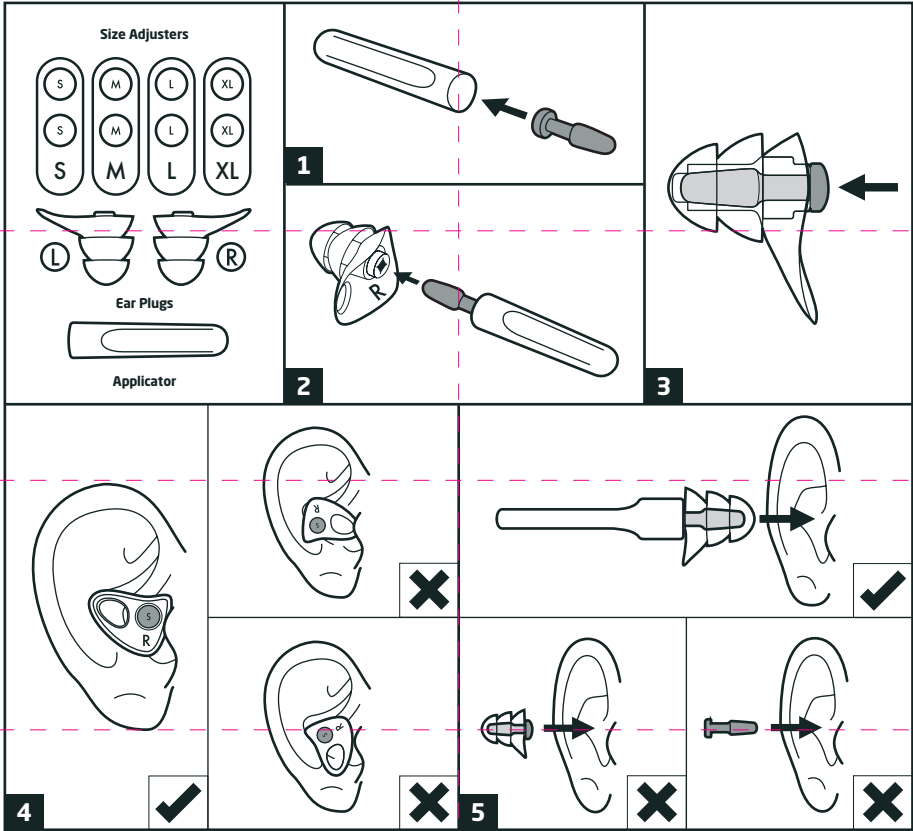
Attenuation Table

Frequency (Hz)	63	125	250	500	1000	2000	4000	8000
Mean Attenuation (dB)	19.0	23.4	23.1	22.3	26.1	28.7	34.4	39.4
Standard Deviation (dB)	5.4	5.6	5.7	2.9	4.0	2.4	4.3	4.9
Assumed Protection (dB)	13.6	17.8	17.4	19.4	22.1	26.3	30.1	34.5

M = 27

M = 22

L = 20

SNR = 26

ENG

Instructions for Use

The Ear Plugs must be fitted correctly for them to work effectively. Inserting anything into your ear must be done slowly and with care. By fitting the different sized Size Adjusters into the Ear Plug silicone ear plugs, one can ensure they fit most ear canals perfectly and this process only needs to be done once.

How To Fit

- Take the Ear Plug for the right ear denoted by the letter "R" printed on the inside of the ear plug.
- Gently push the small Size Adjuster into the Applicator (See Fig. 1) and insert into the Ear Plug (See Fig. 2).
- Ensure each Size Adjuster is pushed as deep into the Ear Plug as possible before fitting them into the ear (See Fig. 3).
- Double checking the orientation is correct (See Fig. 4). Take the Applicator with the Ear Plug attached and carefully insert into the ear canal (a gentle pull of the ear with the opposite hand opens the ear canal making insertion easier). Carefully pull back the Applicator to leave the Ear Plug perfectly fitted inside the ear.
- Repeat steps 1-4 with the left Ear Plug denoted by the letter "L".
- Once fitted, brush your palm over your ears. If you feel the Ear Plug's edges, it is possible the orientation is not correct. Check and re-insert if necessary.

Enlarging the Ear Plug to Improve Fitting

If the Ear Plug feels too loose in the ear canal, the small Size Adjuster may be too small and a larger size may be preferred. Follow the steps below to change the Size Adjuster to your preferred size.

- Pull the small Size Adjusters with fingers from both Ear Plug to remove. Note: Ear Plugs are made from highly durable silicone and are extremely elastic. They will stretch significantly when removing the Size Adjusters.
- Follow steps 1-5 of HOW TO FIT to insert the medium, large or extra large Size Adjuster until the fit is comfortable.

Care Instructions

- These Ear Plug are made from durable silicone and therefore should be hand washed in mild detergent liquid to remove any residue.
- To disinfect the ear plugs wash in a mild disinfectant and wash thoroughly. This product may be adversely affected by certain chemical substances. Further information should be sought from the manufacturer.

Recommendations

- Ensure that the Ear Plugs are fitted, adjusted and maintained in accordance with these instructions.
- Ensure that the Ear Plugs are worn at all times in noisy surroundings.
- For temporary use only.
- In the event they are used daily for 4 months, they must be replaced.

Storage

- The Ear Plugs should be cleaned after use and kept in the carry case supplied.

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